

VAMP Jr. System Inservice Bulletin


Venous Arterial blood Management Protection system

Safe and Accurate Closed Blood Sampling System


Featuring the Edwards TruWave Disposable Pressure Transducer

PRIMING


- 1 Open reservoir plunger to approximately 1/2 cc to facilitate flow of priming solution.

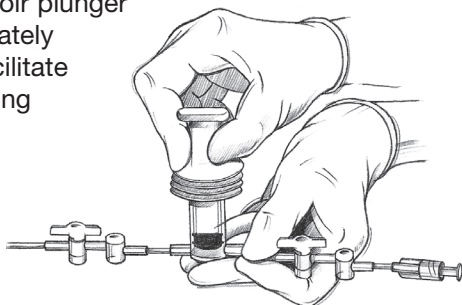


- 2 With the shut-off valve in the open position (parallel to the tubing), hold sampling site above the Edwards VAMP Jr. reservoir at 45° angle.

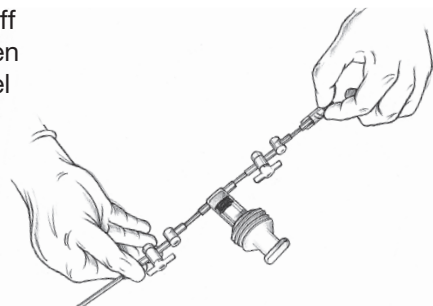

- 3 Provide flow by pulling Snap-Tab of the Edwards TruWave disposable pressure transducer. Slowly deliver priming solution to remove air.

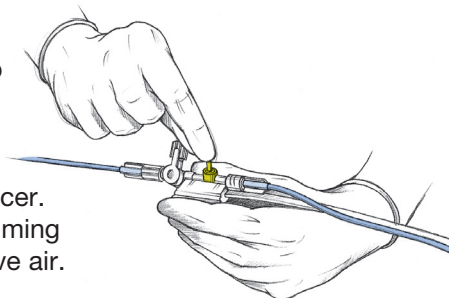


- 4 Close plunger and remove all air from the system. Ensure that VAMP Jr. and catheter connection is a fluid-to-fluid connection.

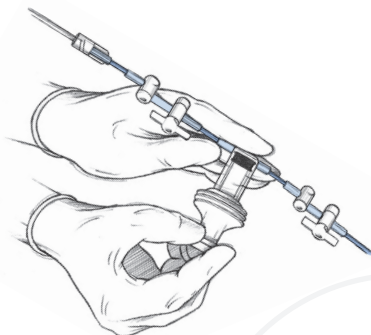


DRAWING THE CLEARING VOLUME

- 5 Close the distal shut off valve by turning the handle perpendicular to the tubing. Smoothly and evenly pull up on the reservoir plunger to draw the required amount of clearing volume consistent with the patient's clinical condition, patient's size, or your hospital protocol. (Recommended pull rate is 1 cc every 10-15 seconds.)

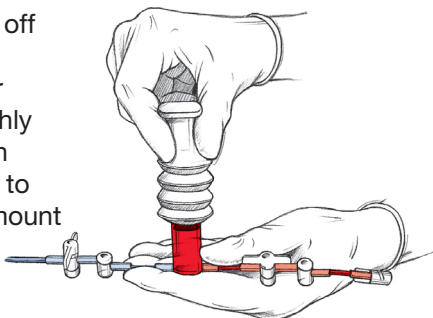


- 6 Close the proximal shut-off valve by turning handle perpendicular to tubing.

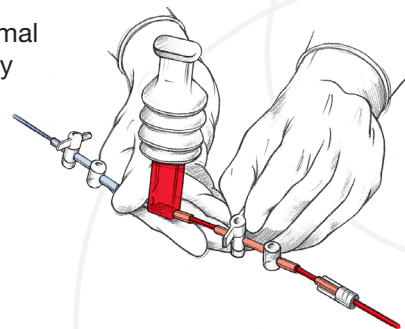


DRAWING BLOOD SAMPLES FROM THE VAMP JR. NEEDLELESS SAMPLING SITE

- 7** Swab proximal sample site. Ensure plunger is depressed to the bottom of the syringe.


Attach cannula. Push cannula into proximal sample site.

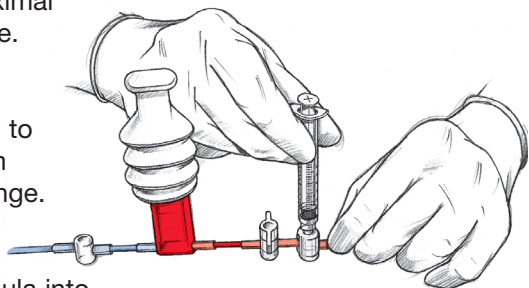


- 8** Slowly draw the blood sample.

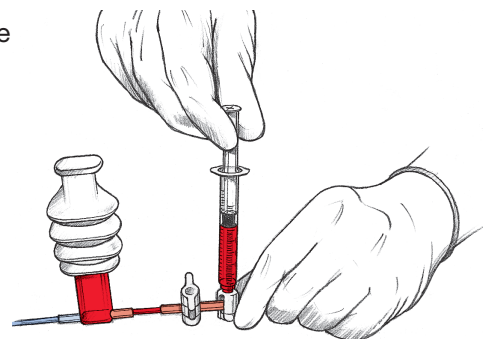


- 9** While holding the cannula, remove the syringe/cannula from the sampling site by pulling straight out.

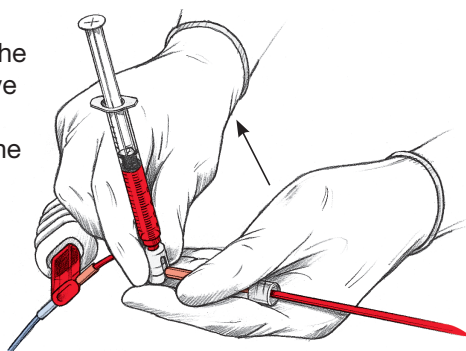


- 11** Slowly, smoothly, and evenly, push down on the plunger until it is fully closed. (Recommended rate is 1 cc every 10-15 seconds.)

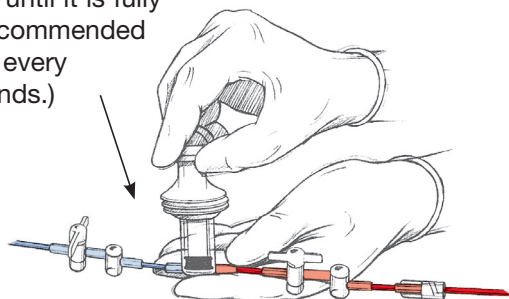

- 10** Open the proximal shut-off valve by turning handle parallel to tubing.

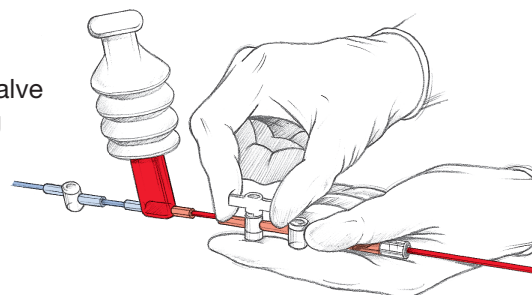

- 12** Fill syringe with flush solution. Attach Interlink NeedleLess cannula to the syringe. Ensure syringe and cannula are free of air bubbles. Swab distal access port. Insert syringe and cannula into the distal access port. Slowly flush line per hospital procedure. Open the distal one-way valve.

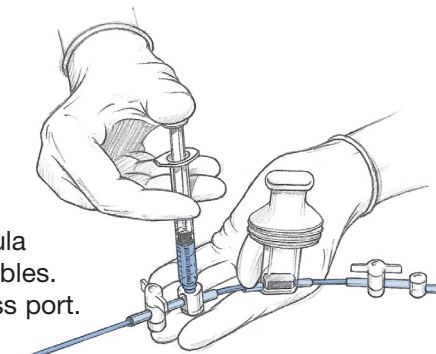

No components of this package or the product it contains are made from natural rubber latex or dry natural rubber.

For professional use. CAUTION: Federal (United States) law restricts this device to sale by or on the order of a physician. See instructions for use for full prescribing information, including indications, contraindications, warnings, precautions and adverse events.

Edwards is a trademark of Edwards Lifesciences Corporation. Edwards Lifesciences, the stylized E logo, TruWave, VAMP and VAMP Jr. are trademarks of Edwards Lifesciences Corporation and are registered in the United States Patent and Trademark Office.

© 2011 Edwards Lifesciences Corporation. All rights reserved. AR05992

Edwards Lifesciences

Irvine, USA | Nyon, Switzerland | Tokyo, Japan | Singapore, Singapore | São Paulo, Brazil

edwards.com

